


MAb : Mouse monoclonal IgG1
(anti-human IgG4)

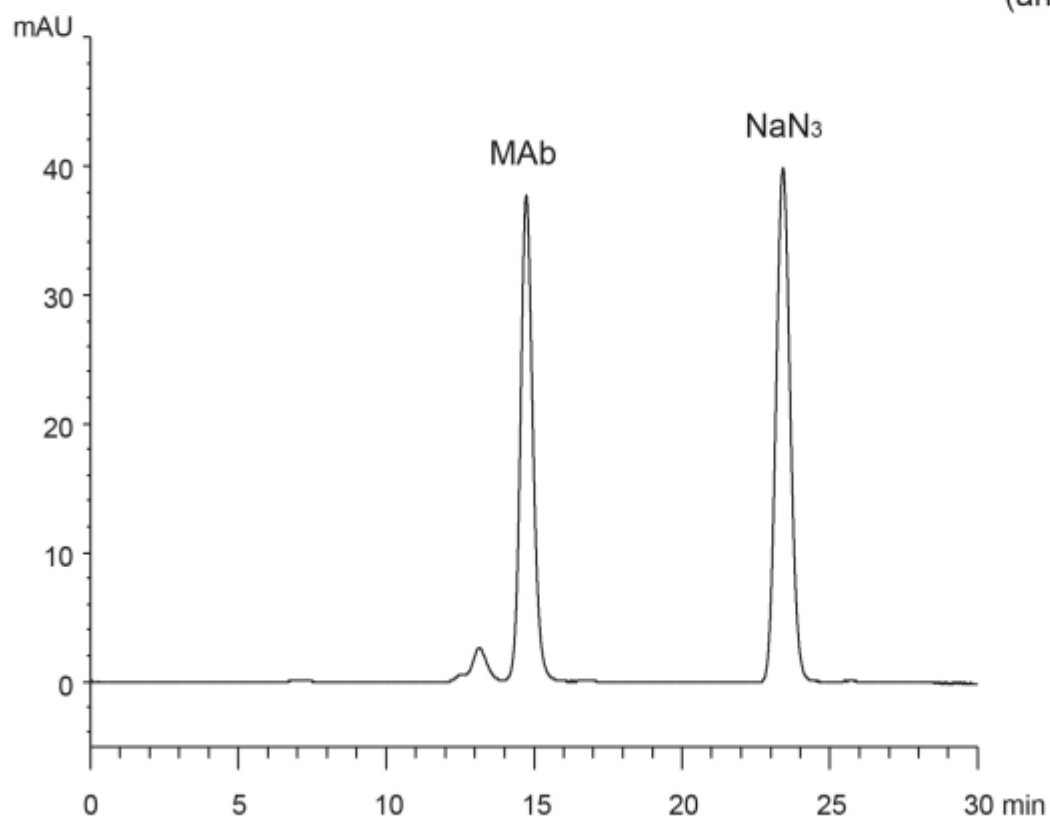

Column : YMC-Pack Diol-200
300 X 4.6 mm I.D.

Eluent : 0.1 M KH_2PO_4 - K_2HPO_4 (pH 7.0)

Flow rate : 0.17 mL/min

Temperature : ambient (25°C)

Detection : UV at 220 nm

Injection : 10 μL

Sample : a commercially available mouse monoclonal IgG1 (0.05 mg/mL)
(purified by DEAE chromatography, containing NaN_3)